

PHERO-sub 0710-E




- Injection moulded body
- Easy lid removal
- Available with 7 x 7 cm geltray
- Available with 7 x 10 cm geltray
- Available with both geltrays

The PHERO-sub 0710-E is the smallest unit in the range, designed for low medium numbers of samples. The smallest gel size maximises run economy but does not compromise versatility as two tray options are available 7 x 7 cm and 7 x 10 cm, and combs ranging from preparative up to 16 samples. Simply by altering the gel tray or comb, this compact unit is capable of resolving up to 64 different samples, prepping 1 ml of sample or separating sample bands over a distance of 9 cm.

Technical specifications

<i>Gel dimensions (W x L)</i>	<i>7 x 7 cm 7 x 10 cm</i>
<i>Unit dimensions</i>	<i>21 x 9 x 9 cm</i>
<i>Max sample capacity</i>	<i>7 x 7 cm – 32 samples 7 x 10 cm – 64 samples</i>
<i>Buffer volume</i>	<i>225 ml</i>
<i>Combs</i>	<i>1, 2, 4, 8(MC), 8, 10, 12(MC), 16</i>
<i>Thicknesses</i>	<i>0.75, 1, 1.5, 2 mm</i>

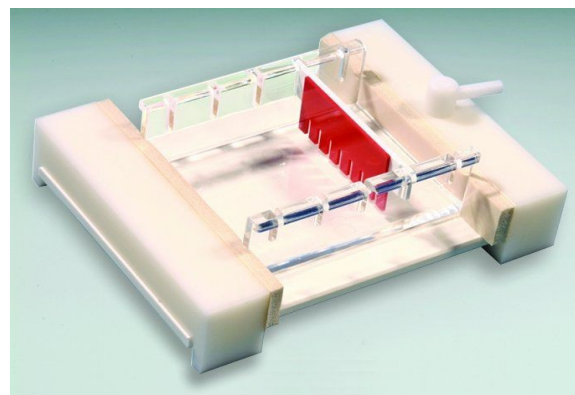

Extended sub flexi caster