




D I E K R E A T I V E A N A L Y S E N T E C H N I K

Model 1000.

The intelligent solution for your electrophoresis laboratory.

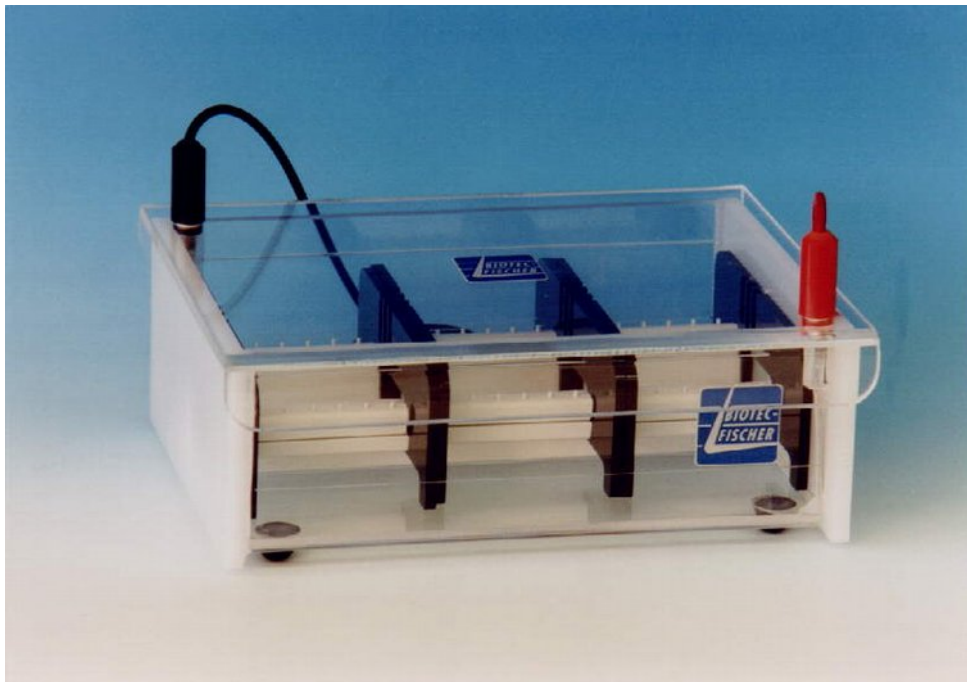

Electrophoresis chamber model 1000 with strip bridge for the PRA application system

- Universal electrophoresis chamber for serum protein, lipoprotein, Hb-A1, haemoglobin electrophoresis and electrophoresis of concentrated urine and cerebrospinal fluid
- For simultaneous separation of up to 24 samples
- Extremely safe, automatic cut-out when the chamber is opened
- Extremely economical with low purchase price and low buffer volume
- Very easy to operate and clean
- Lightweight with minimal space requirements
- Suitable for a variety of application systems

■ **Universal.**

The electrophoresis chamber model 1000 is a universal system for serum protein, lipoprotein, Hb-A1 and haemoglobin electrophoresis as well as electrophoresis of concentrated urine and CSF. Up to 24 samples can be electrophoretically separated simultaneously. It can accommodate all the strip bridges for commercially available application systems (e.g. PRA and CEL, Bender & Hobein, Boskamp, Paragon Im-munfixation, etc.).

■ **Safe.**

As a safety precaution, lifting the lid immediately stops the supply of electric current. This guarantees a high degree of safety for the user. The two connecting cables are permanently fitted to the cover of

the chamber, ensuring that they are never mixed up and the user cannot touch live parts.

■ **Robust.**

The buffer tank and the safety cover of the model 1000 are made from strong, hard-wearing plexiglass. This makes the unit very robust and damage-resistant.

■ **Compact.**

The small-sized, lightweight model 1000 is very easy to operate and will take up very little space in your laboratory.

■ **Economical.**

At a very reasonable purchase price and requiring very little buffer, the model 1000 is extremely economical to use.

■ **Standard delivery details.**

The standard delivery comprises the buffer tank with safety cover and a set of cables. The application system and power supply have to be ordered separately.

■ **Accessories.**

A wide variety of accessories can enhance operating comfort and increase the application options for your unit.

(Specifications may be subject to change)

Ordering information:

<u>Code</u>	<u>Product</u>	<u>Quantity</u>
321001	Electrophoresis chamber model 1000 , standard delivery (complete with safety cover and fitted connecting cables)	1 set